**seniorfon®**

Vivanco GmbH  
Ewige Weide 15  
D-22926 Ahrensburg/Germany  
Fax +49(0)41 02 / 231 - 160  
www.babyfon.com









055

- D** Bedienungsanleitung
- GB** Instruction manual
- F** Mode d`emploi
- E** Instrucciones de uso
- I** Istruzioni per l`uso
- NL** Gebruiksaanwijzing
- SF** N käyttöohjet
- N** Betjeningsveiledning

**BM 770**

<b>D</b>	Bedienungsanleitung . . . . .	4
<b>GB</b>	Instruction manual . . . . .	10
<b>F</b>	Mode d`emploi . . . . .	16
<b>E</b>	Instrucciones de uso . . . . .	22
<b>I</b>	Istruzioni per l`uso . . . . .	28
<b>NL</b>	Gebruiksaanwijzing . . . . .	34
<b>SF</b>	N käyttöohjet . . . . .	40
<b>N</b>	Betjeningsveiledning . . . . .	46

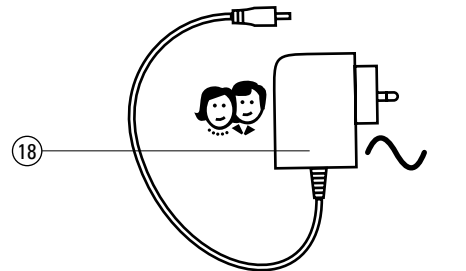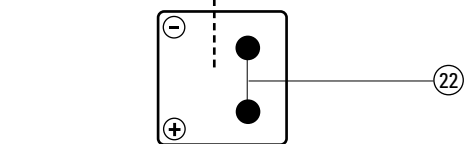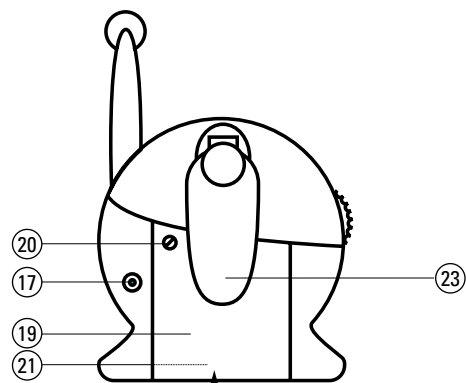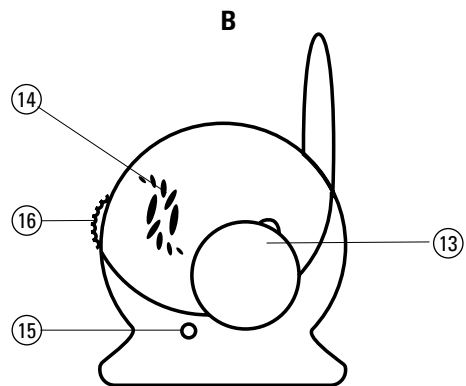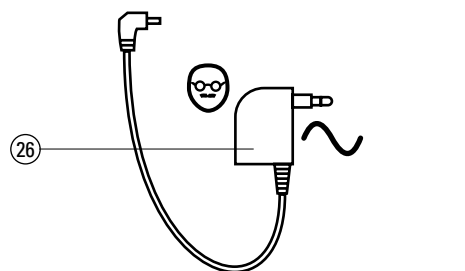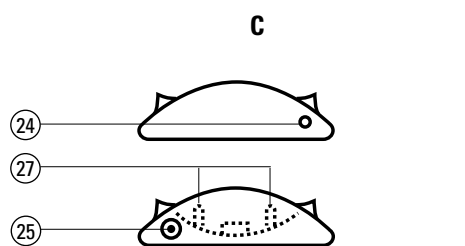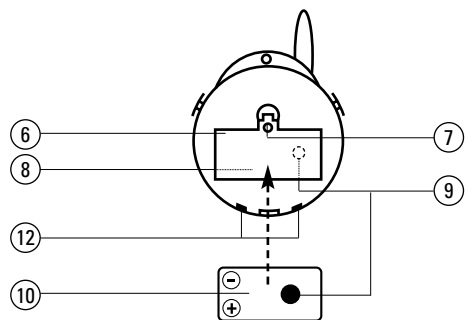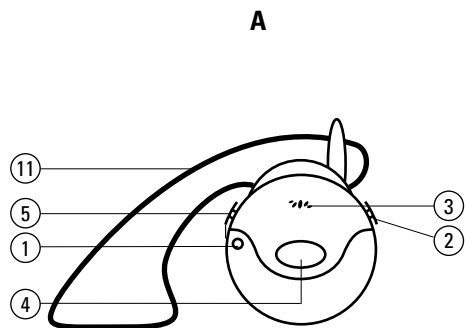**Seniorfon BM 770**

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**BM 770**

Babyfon® is one of the leading brands for wireless baby monitors in Europe. We have devoted ourselves for many years to the development and production of reliable, high-quality products. With the Seniorfon BM 770 we are rounding off our range with a device for the elderly. The Seniorfon BM 770 is fitted with a manual alarm function tailored specially to seniors. It also provides all the familiar functions used in baby supervision.

**A Senior device**

- 1 On/Off/Channel selector switch
- 2 Power-on lamp
- 3 Microphone
- 4 Alarm button
- 5 Sensitivity setting
- 6 Battery compartment lid
- 7 Screw for battery compartment lid
- 8 Battery compartment
- 9 Charge contacts for battery pack
- 10 Battery pack
- 11 Neck strap
- 12 Charge contacts for charger

**B Reception device**

- 13 On/Off/Channel selector switch
- 14 Speaker
- 15 Power-on lamp
- 16 Volume control
- 17 AC adapter connector
- 18 AC adapter
- 19 Battery compartment lid
- 20 Screw for battery compartment lid
- 21 Battery compartment
- 22 Charge contacts for battery pack
- 23 Belt clip


**C Charger for senior device**

- 24 Power-on indicator
- 25 AC adapter connector
- 26 AC adapter
- 27 Charge contacts for senior device


**Important information**

Please read these instructions carefully. This Seniorfon is only intended as an aid. It is not a substitute for responsible, correct supervision and should therefore not be used as such.


- Never use the Seniorfon near water or in damp places.
- Avoid exposure to extreme heat or strong light.
- Never open the housings or AC adapters (danger of electric shock!).
- Make sure that the local mains voltage corresponds to the voltage specified on the AC adapter.
- Avoid touching the battery and charging contacts. Do not touch the contacts with metal objects.
- Use a damp cloth to clean the Seniorfon.

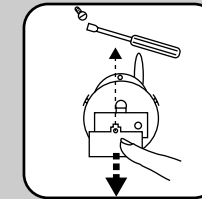**Preparing your Seniorfon****Senior device**


1. Open the battery compartment by loosening the screw on the battery compartment lid.

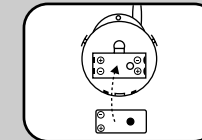
2. Place the battery pack supplied in the battery compartment, observing the correct polarity. Make sure that the charge contact in the battery compartment is established (see A9 on senior device).

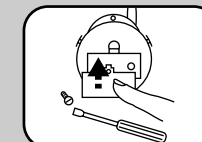
3. Slide the battery compartment lid back on to the battery compartment and tighten the screw.

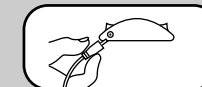
4. Connect the AC adapter to the connector of the charger for the senior device. Place the senior device in the charger so that the charge contacts of the senior device and the charger are connected.

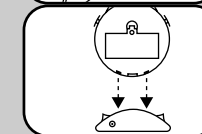5. The battery pack is now charged automatically. Charging is complete after approx. 15 hours.

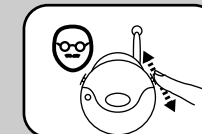6. Switch on the senior device by selecting channel A or B.

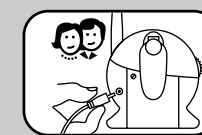7. Take the reception device and move away from the senior device. If the units are placed too close together, there may be a high-pitched whistling sound.



8. Connect the second AC adapter to the connector of the reception device.

9. Switch on the reception device and select the same channel as on the senior device.

10. Your Seniorfon is now ready for use. Transmission begins automatically when the senior device is activated by a noise.
11. Make a noise to test whether transmission is ensured at the location you have chosen.

### Using your Seniorfon

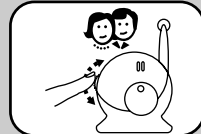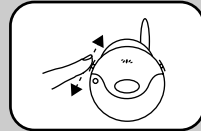1. The senior device can be operated from a fixed position in the charge cradle. It has all the functions you are familiar with from a normal radio frequency Babyfon. As soon as the senior device is activated by a noise, it transmits the signals to the reception device.
2. The senior device can also be operated as a mobile unit, i.e. persons requiring supervision can wear the senior device using the neck strap.
3. Make sure that both devices are switched on and that the same channel is selected.
4. Set the response sensitivity of the microphone on the senior device to the level you require. To begin with, we recommend a middle position that can later be adjusted as required.
5. Set the volume on the senior device to the level required. If the senior device is activated by a noise, transmission starts automatically and the power-on lamp of the reception device changes from green to red for the duration of the transmission.

### Explanation of functions

#### Alarm button on the senior device

The senior device is fitted with an alarm button, enabling the person wearing it to activate a manual alarm at the touch of a button. An alarm tone then sounds for approx. 5 seconds on the reception device.

#### Reception monitoring

There is regular radio contact between the senior device and the reception device. A test is made as to whether there is still sufficient contact to the senior device. If radio contact breaks down completely, for example because the senior device is outside the range, or because the batteries in the senior device are empty, the reception device emits an alarm tone after a short time. In addition, the power-on lamp flashes. The alarm tone and flashing stop automatically soon as radio contact is resumed.

#### Digital pilot tone / preventing interference

Besides sounds, the senior device also transmits an inaudible signal tone. The reception device only begins reception when this signal tone is received. This practically rules out the reception of external wireless systems. In rare cases, however, interference from other wireless devices in the vicinity may occur. In the event of such interference or poor reception, switch both devices to the other channel.

#### 250 m range


In ideal conditions, the range of your Seniorfon outdoors is approx. 250 metres. Obstacles such as walls, ceilings or large objects affect the range. The more obstacles there are between the senior and reception devices, the smaller the range becomes.

#### Battery pack operation of the reception device (battery pack not included)

The supplied battery pack can only be used with the senior device. You can order a battery pack (Babyfon® model number EAB 4) for the reception device from a specialist dealer. This battery pack can then also be charged using the charge function integrated in the reception device.

#### Using other battery packs or batteries

The reception and senior devices can also be powered with standard AAA battery packs or non-rechargeable 1.5V AAA batteries. However, standard AAA battery packs cannot be charged using the devices.
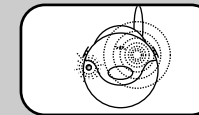
#### Battery alarm

A signal tone and the flashing power-on lamp warn you as soon as the battery pack or batteries in the senior device are empty. You should then insert new batteries or recharge the battery pack.

#### Disposal of batteries

Consider the environment when disposing of the batteries, which must not be thrown away with domestic rubbish. Do not use old and new batteries or different battery types simultaneously. Remove batteries from the device if it is not to be used for a long time.

### Specifications

No. of channels:	2	
Frequency	864.8 – 865.0 MHz	
Power supply:	Receiver:	AC adapter plug 9 Volt/ 200mA or battery pack (4x AAA batteries) (not included, available from specialist dealers under Babyfon model number EAB 4)
	Transmitter:	AC adapter plug 9 Volt/ 200mA or battery pack EAB 2 (2x AAA batteries) (included)
Transmitting power:	max. 10 mW	
Operating temperature:	10 to 40° C	

Subject to alteration.

Babyfon® is a registered trademark of Vivanco AG, Ahrensburg

### Trouble-shooting

If your Seniorfon BM 770 should not function properly, please proceed as follows:

1. Using the following table of malfunctions, try to find and correct the fault.
2. Read the corresponding section in these instructions carefully.
3. Look on the Internet at [www.babyfon.de](http://www.babyfon.de) or [www.vivanco.com](http://www.vivanco.com) to see whether the fault and its correction are described.
4. Ask your dealer.
5. Call the telephone hotline 01805 / 40 49 10 (€0.12/min.) for Germany. You will find contacts for Europe at [www.vivanco.com](http://www.vivanco.com).

<b>Problem</b>	<b>Possible cause</b>	<b>Solution</b>
No transmission / no sound	Different channels set on both devices	Set both devices to same channel (A1/B13)
	Transmitter does not switch automatically to transmit	Increase response sensitivity (A5)
	Devices are outside range	Reduce distance between devices
	Transmitter is not supervised person or too far away	Reduce distance of transmitter from supervised person
	Batteries are empty	If batteries are used, replace them
Interference or another voice can be heard during transmission	Unfavourable choice of channel	Change to another channel (A1/B13)
Loud beeping and whistling noise	Feedback due to insufficient distance between devices	Increase distance between devices or reduce volume on receiver (B16)
Battery packs no longer being charged	AC adapter plug is not connected to power supply	Check connection
	Plug of AC adapter not connected properly to device	Insert plug properly into connector socket (B17/C25)
	Senior device not positioned correctly on charge cradle	Position senior device correctly until power-on indicator glows
	Battery packs have come to the end of their service life	Replace battery packs
Seniorfon transmits every tiny sound	Response sensitivity set too high	Reduce response sensitivity (A5)


**Licence**

**CE 0700**

This Vivanco FM product conforms with the European R&TTE directive for FM systems in harmonised frequency ranges. These systems can be sold and operated in the EU, as well as in CH and N.  
*Operation is not permitted in all countries not named.*  
 The declaration of conformity is included at the end of this manual.

**Warranty**

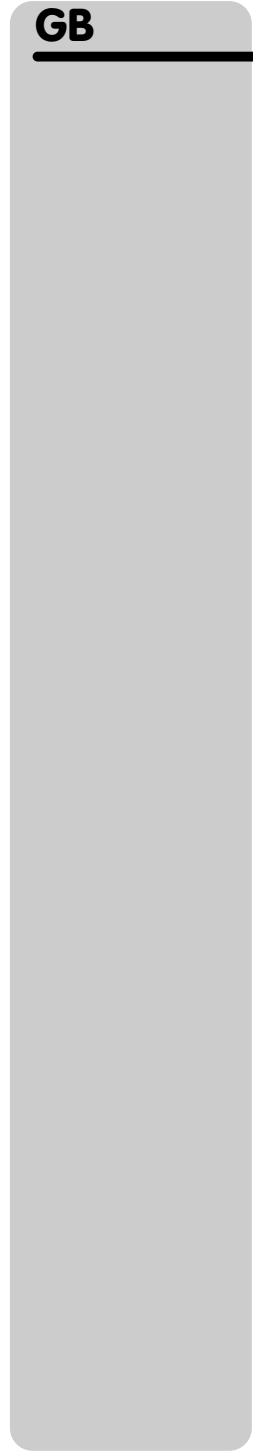The warranty period is 12 months from the date of purchase. Please contact your specialist dealer or  
 Vivanco Service  
 Ewige Weide 17,  
 D-22926 Ahrensburg/Germny  
 A copy of the purchase receipt must be included.

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Babyfon® er et registrert varemerke som tilhører Vivanco AG, Ahrensburg